



EDOF IOLs

Active

Active SL (SeeLens MF)

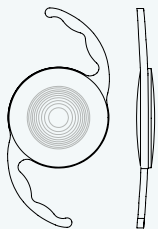

* Active BN EZ (BunnyLens MF EASY)

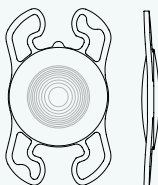


* Available in preloaded configuration

	Active SL (SeeLens MF)	Active BN EZ (BunnyLens MF EASY) *
Platform	C-Loop	4 - Loop
Overall Diameter	13 mm	11 mm ($\geq 10D$) 11.5mm ($< 10D$)
Placement	Capsular Bag	
Optic Diameter	6 mm	
Power Range	0 to +30 (0.5D increments) +30 to +35 (1D increments)	
Addition	Add power: 2.25	
Cylinder Range	-	-
Optic Design	Aspheric Multifocal Diffractive Apodized	
360° Continuous Square Edge	Yes	
Haptic Angulation	5°	
Material	Hydrophilic acrylic with bonded UV absorber and violet light filter	
Refractive Index	1.46 (hydrated @ 35°C)	
A-constant (SRK/T) for Optical or Immersion US Biometry	118.6	118.5
A-constant (SRK/T) for Contact US Biometry	118.26	118.16
Sterilization	Steam	

Provide Your Patient With a Real EDof Solution

The Active lens employs a diffractive pattern that allows your patient to be spectacle-free for both far and intermediate visual ranges.


The Active Defocus Curve

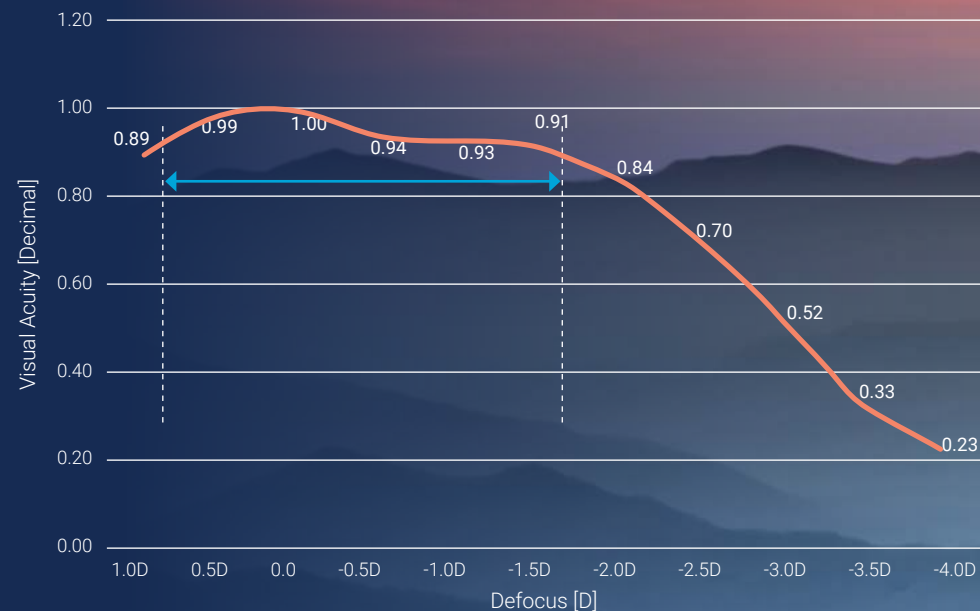

Study by Joon Hyun Kim MD, PhD; ESCRS 2018

Sharpest Square Edge

The novel manufacturing technology enables to achieve exceptionally sharp 360° square edge, proven to be effective against PCO along with a wide-angle contact with the capsular bag.

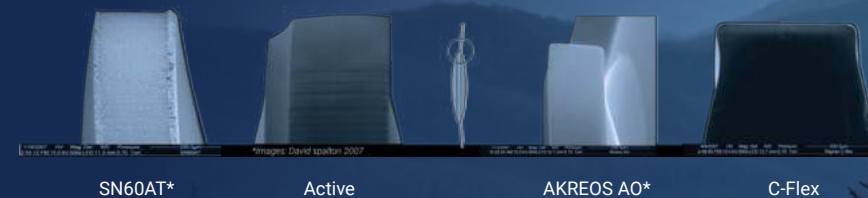

* Seelens MF / Bunnylens MF Easy +2.25